

stryker®

SILVERGlide® PRO SERIES

Neuro Spine ENT



OPTIMIZED
FOR VISIBILITY

SMIS

SILVERGlide® PRO SERIES

The MIS Focus Section

- Ultra-thin tines allow for smaller incisions and the ability to gain deeper access in delicate applications
- Minimizing the forceps design and reducing weight provides better balance and reduced fatigue during lengthy procedures



Efficiency

- Reduced scissoring feature allows
 - faster, continual coagulation of targeted tissue
 - better grasp and manipulative control
- Thinner tines, longer length and a slimmer handle provide better visibility for MIS cases




ProSeries Features

- **Handle:** Light-weight, easy to grip and manipulate
- **Guide-stop:** Constant parallel alignment of the hardened tines increases control and grasping ability when manipulating tissue.
- **Tines:** Ultra-thin hardened tines allow better access and improved line of sight in virtually any surgical application regardless of incision size or depth of the wound site.
- **Tips:** Silver plating assists in dissipating heat quickly, reducing collateral heating and hot spots that cause sticking.




Thermal Spread

The thermal image illustrates that the silver plated tips reduce thermal spread to adjacent tissue. The energy is focused on target tissue in between the tips. In the image, yellow denotes warmer temperatures while blue denotes cooler temperatures. The targeted tissue in between the tips is coagulated at higher temperatures, while the thermal spread in adjacent tissue remains low and cooler than temperatures in between the tips. The silver plating efficiently delivers and dissipates heat quickly preserving the healthy margins of the adjacent tissue as well as reducing hot spots and sticking.

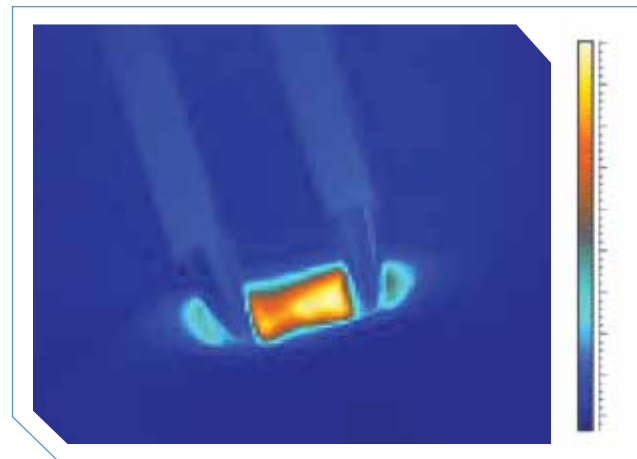

Productivity

Silver plated tips reduce hot spots that cause sticking, reducing the need to continually stop and wipe off eschar saving time and frustration throughout the case.

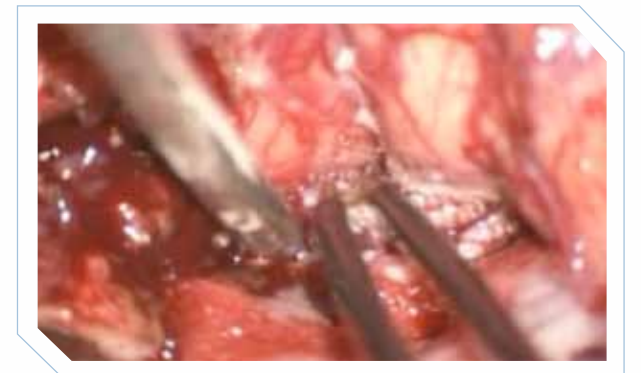

To place an order:
 Call Stryker Customer Service at 800-253-3210
 or contact your Stryker NSE sales representative



ProSeries Order Information and Part Numbers

Part Number	Description
● 14.5cm Bayonet 0.4mm Tip PRO-Series	6730-145-004
● 14.5cm Bayonet 0.7mm Tip PRO-Series	6730-145-007
● 14.5cm Bayonet 1.0mm Tip PRO-Series	6730-145-010
● 14.5cm Bayonet 1.2mm Tip PRO-Series	6730-145-012
● 14.5cm Bayonet 1.5mm Tip PRO-Series	6730-145-015
● 16cm Bayonet 0.4mm Tip PRO-Series	6730-160-004
● 16cm Bayonet 0.7mm Tip PRO-Series	6730-160-007
● 16cm Bayonet 1.0mm Tip PRO-Series	6730-160-010
● 16cm Bayonet 1.2mm Tip PRO-Series	6730-160-012
● 16cm Bayonet 0.4mm, 45° Tip up PRO-Series	6730-160-504
● 19cm Bayonet 0.4mm Tip PRO-Series	6730-190-004
● 19cm Bayonet 0.7mm Tip PRO-Series	6730-190-007
● 19cm Bayonet 1.0mm Tip PRO-Series	6730-190-010
● 19cm Bayonet 1.2mm Tip PRO-Series	6730-190-012
● 19cm Bayonet 1.5mm Tip PRO-Series	6730-190-015
● 19cm Bayonet 2.0mm Tip PRO-Series	6730-190-020
● 19cm Bayonet 0.4mm, 45° Tip up PRO-Series	6730-190-504
● 19cm Bayonet 1.2mm, 45° Tip up PRO-Series	6730-190-512
● 19cm Bayonet 1.2mm, 45° Tip down PRO-Series	6730-190-612
● 22cm Bayonet 0.4mm Tip PRO-Series	6730-220-004
● 22cm Bayonet 0.7mm Tip PRO-Series	6730-220-007
● 22cm Bayonet 1.2mm Tip PRO-Series	6730-220-012
● 22cm Bayonet 1.5mm Tip PRO-Series	6730-220-015
● 22cm Bayonet 2.0mm Tip PRO-Series	6730-220-020
● 22cm Bayonet 0.7mm, 15° Tip up PRO-Series	6730-220-107
● 22cm Bayonet 0.7mm, 30° Tip up PRO-Series	6730-220-307
● 22cm Bayonet 0.7mm, 45° Tip up PRO-Series	6730-220-507
● 22cm Bayonet 1.2mm, 45° Tip up PRO-Series	6730-220-512
● 22cm Bayonet 0.7mm, 45° Tip down PRO-Series	6730-220-607
● 22cm Bayonet 1.2mm, 45° Tip down PRO-Series	6730-220-612
● 24cm Bayonet 0.4mm Tip PRO-Series	6730-240-004
● 24cm Bayonet 0.7mm Tip PRO-Series	6730-240-007
● 24cm Bayonet 1.2mm Tip PRO-Series	6730-240-012
● 24cm Bayonet 1.5mm Tip PRO-Series	6730-240-015
● 24cm Bayonet 2.0mm Tip PRO-Series	6730-240-020
● 24cm Bayonet 0.7mm, 45° Tip up PRO-Series	6730-240-507
● 24cm Bayonet 0.7mm, 45° Tip down PRO-Series	6730-240-607
● 24cm Bayonet 0.7mm, 80° Tip up PRO-Series	6730-240-807
● 24cm Bayonet 0.7mm, 80° Tip down PRO-Series	6730-240-907
● 28cm Bayonet 1.2mm Tip PRO-Series	6730-280-012
● 28cm Bayonet 1.5mm Tip PRO-Series	6730-280-015
● 19cm Straight 0.7mm Tip PRO-Series	6735-190-007
● 19cm Straight 1.0mm Tip PRO-Series	6735-190-010
● 19cm Straight 1.0mm, 30° Tip up PRO-Series	6735-190-310
● 19cm Straight 2.0mm, 30° Tip up PRO-Series	6735-190-320

Reconstructive

- Hips
- Knees
- Trauma & Extremities
- Joint Preservation
- Orthobiologics

Medical & Surgical

- Power Tools & Surgical Accessories
- Image Guided Navigation
- Endoscopy & Arthroscopy
- Integrated Communications
- Beds, Stretchers & EMS
- Sustainability Solutions

Neurotechnology & Spine

- Cranio-maxillofacial
- Interventional Spine
- Neurosurgical, Spine & ENT
- Neurovascular
- Spinal Implants

Stryker® Neuro Spine ENT
 4100 East Milham Avenue
 Kalamazoo, MI 49001 USA
 t: 269 323 7700 f: 269 353 3026
 toll free: 800 253 3210
 www.stryker.com

Products referenced with TM designation are trademarks of Stryker. Products referenced with ® designation are registered trademarks of Stryker.

Literature Number: 9100-001-775 Rev. None
 UnDe/P.S.
 Copyright © 2012 Stryker
 Printed in USA

